

# SYST'AM® P101M / POLYPLOT

A MOULDED WAFFLE MATTRESS MADE OF HR FOAM WITH HIGH DENSITY

## MATERIALS

40  
kg/m<sup>3</sup>



HR foam



PUR POLYFILM  
fitted bed sheet



POLYMAILLE®  
integral cover



POLYMAILLE®**HD**  
integral cover



NF EN ISO 597 - 1 & 2

### Foam maintenance:



(Do not immerse in water)

### Cleaning POLYFILM cover:



### Cleaning POLYMAILLE® cover:




### Cleaning POLYMAILLE® HD cover:



## INDICATIONS

### Risk

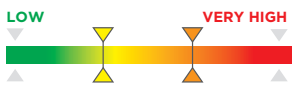

P101M / POLYPLOT Mattress can be used, when:

- 1- Confinement in bed is necessary resulting from diseases or handicap and there is at the same time an increased risk of developing pressure ulcer according evaluation scales (Braden, Norton, Waterlow, ...).
- 2- pressure ulcer (Stages 1 & 2 according EPUAP) has developed.



### DELIVERED COMPRESSED AND ROLLED

→ Gain in storage volume 2.5 time less volume.



14 cm

### ROUNDED CORNERS

→ Easier installation of supports in bed base.



### PERFORMANCE AND COMFORT OF THE SYST'AM® HIGH RESILIENCE FOAM

- Because it is moulded and of high density (40 kg/m<sup>3</sup>), it ensures durable flexibility, comfort and a soft touch.
- It offers as well a better resistance against external wear and tear, for a longer lifespan of the mattress.

### MATTRESS MOULDED IN ONE SECTION

- Guarantee for a perfect holding in the cover.
- Compatible with standard and multi-sections bed bases.
- Gain of time and easier installation of the supports.



### LARGE SURFACE OF PLOTS

→ Allowing to maximise the contact surface between the patient and the support thus to decrease the Max. pressures on area at risk.



### INDEPENDENCE & MOBILITY OF THE BLOCKS

→ Have been carefully researched to offer a significant reduction in the effects of chafing.



### AERATION CHANNELS BETWEEN BLOCKS

→ Help to reduce the effects of maceration through improved mattress ventilation.

## AVAILABLE COVERS






POLYFILM FITTED BED SHEET  
FOR SINGLE PATIENT USE

- A breathable and hyper-extensible hypoallergenic polyurethane film.
- Unfailing adaptability: Unlike PVC, this material is extremely strong and offers lasting resistance against wear and tear.
- POLYFILM is totally impermeable to liquids.



POLYMAILLE® & POLYMAILLE® **HD**  
INTEGRAL COVER FOR  
MULTI-PATIENT USE

- Hypoallergenic polyurethane enduction (anti-bacterial, anti-germ, anti-fungal).
- POLYMAILLE® is bi-extensible, to reduce tear and shear forces.
- Zip over 3 sides. Non-slip base (POLYMAILLE® only).
- POLYMAILLE®**HD**: welded seams

DESIGNATION	ITEM CODE	 (cm / inches)	 (cm / inches)	 min - MAX (kg / lbs)
P101M / POLYPLOT with POLYFILM fitted bed sheet	P101MPLPFDH	87 x 198 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286
P101M / POLYPLOT with POLYMAILLE® integral cover	P102MPLPMHI	87 x 198 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286
	P102MPLPMHI190	87 x 187 / 34 x 74	90 x 190 / 35 x 75	35 - 130 / 77 - 286
	P102MPLPMHI220	87 x 220 / 34 x 87	90 x 220 / 35 x 87	35 - 150 / 77 - 330
	P102MPLPMHI100	96 x 198 / 38 x 78	100 x 200 / 39 x 79	35 - 150 / 77 - 330
	P102MPLPMHI110	103 x 198 / 44 x 78	110 x 200 / 43 x 79	135 - 270 / 297 - 595
	P102MPLPMHI120	113 x 198 / 47 x 78	120 x 200 / 47 x 79	135 - 270 / 297 - 595
P101M / POLYPLOT with POLYMAILLE® <b>HD</b> integral cover	P102MPLPMHI140	133 x 199 / 52 x 78	140 x 200 / 55 x 79	135 - 270 / 297 - 595
	P102MPLPMHF	87 x 199 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286
	P102MPLPMHF190	87 x 187 / 34 x 74	90 x 190 / 35 x 75	35 - 130 / 77 - 286

All mattresses are delivered compressed. Should you want  
**NOT COMPRESSED MATTRESSES**, please add NC at the end of the article-N°



WARRANTY

2  
YEARS