

SYST'AM® P181M / POLYMULTI® VISCO

A MOULDED MATTRESS IN ONE SECTION MADE OF HR FOAM WITH HIGH DENSITY AND 2 VISCOELASTIC FOAM PADS

MATERIALS



Foam maintenance:




Cleaning POLYFILM cover:



Cleaning POLYMAILLE® cover:



INDICATIONS

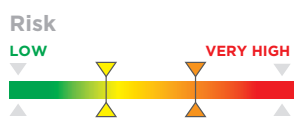

P181M / POLYMULTI® VISCO Mattress can be used, when:

- 1- Confinement in bed is necessary resulting from diseases or handicap and there is at the same time an increased risk of developing pressure ulcer according evaluation scales (Braden, Norton, Waterlow, ...).
- 2- pressure ulcer (Stages 1 & 2 according EPUAP) has developed.



VISCOELASTIC FOAM PADS

- Precise moulding of the body and increase of the body surface in contact with the mattress,
- Better pressures distribution: reduction of transcutaneous pressures on areas at high risk and facilitated blood circulation,
- Improved comfort and stability of the patient,
- Multibearing surface.
- 4 cm of viscoelastic foam pad under the sacral area and the heels area.

AVAILABLE VERSIONS





SYST'AM® P181M / POLYMULTI® VISCO



PERFORMANCE AND COMFORT OF THE SYST'AM® HIGH RESILIENCE FOAM

- Because it is moulded and of high density (40 kg/m³), it ensures durable flexibility, comfort and a soft touch.
- It offers as well a better resistance against external wear and tear, for a longer lifespan of the mattress.

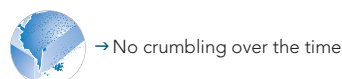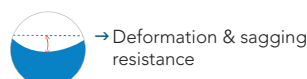

MATTRESS MOULDED IN ONE SECTION

- Guarantee for a perfect holding in the cover.
- Compatible with standard and multi-sections mattress bases.
- Gain of time and easier installation of the supports.

AVAILABLE COVERS



- A breathable and hyper-extensible hypoallergenic polyurethane film.
- Unfailing adaptability: Unlike PVC, this material is extremely strong and offers lasting resistance against wear and tear.
- POLYFILM is totally impermeable to liquids.



- Hypoallergenic polyurethane enduction (anti-bacterial, anti-germ, anti-fungal).
- POLYMAILLE® is bi-extensible, to reduce tear and shear forces.
- Zip over 3 sides. Non-slip base.

DESIGNATION	ITEM CODE	(cm / inches)	(cm / inches)	min - MAX (kg / lbs)
P181M / POLYMULTI® VISCO with POLYFILM bed sheet	P181MPPFDH	87 x 198 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286
P181M / POLYMULTI® VISCO with POLYMAILLE® cover	P181MPPMH	87 x 198 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286

All mattresses are delivered compressed. Should you want **NOT COMPRESSED MATTRESSES**, please add NC at the end of the article-N°

